


SPINEBOARDS

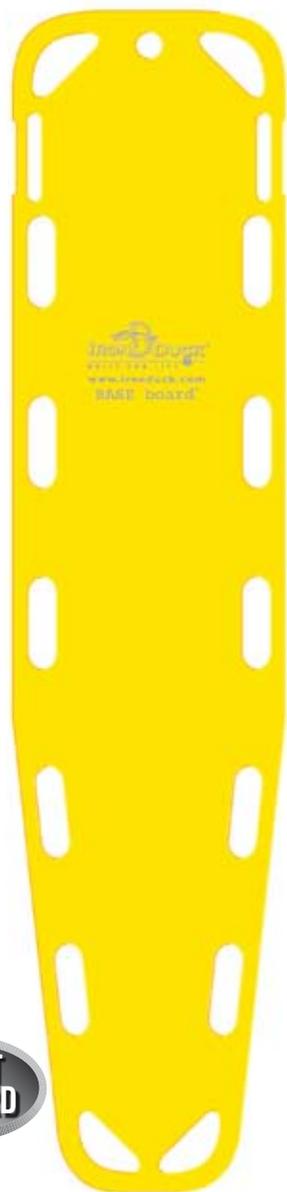

FLAT BOARD

• **Size:** 72" x 16" x 1.75"

• **Weight:** 16 lbs.

• **Capacity:** 500 lb.

• **Float Capacity:** 250 lbs.

Part# 35850

Base Board

- Offers the leading strength and safety in its weight class.
- Rated for the lowest amount of deflection in the 400 to 500 lbs. weight class.
- X-Ray and CT scan translucent, MRI Scan compatible.
- Flat surface patient area won't compromise CPR efforts or aggravate cervical spine injuries.
- Made of a roto-molded HDPE (high-density polyethylene) shell over non-toxic, eco-friendly foam for durability and easy cleaning.
- 14 oversized hand holds can accommodate gloved hands.
- 12 optional offset composite pins can be located in the lower third of the handholds to keep pins out of the way while lifting patients.

COLORS:

- Yellow • Red • Orange • Lime Green • Blue
- Olive Drab • Black

OPTIONS:

- Available with or without 12 board pins
- Custom graphics available. Additional fees apply.



BUILT IN AMERICA



413-377-0030
www.ironduck.com

